

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)


(51) International Patent Classification ⁶ : H01J 49/40, 49/42	A1	(11) International Publication Number: WO 99/30350 (43) International Publication Date: 17 June 1999 (17.06.99)
(21) International Application Number: PCT/CA: (22) International Filing Date: 3 December 1998 (6) (30) Priority Data: 60/067,570 5 December 1997 (05.12.97) (71) Applicant (for all designated States except US): UNIV. OF BRITISH COLUMBIA (CA/CA); 2036 Me. Vancouver, British Columbia V6T 121-(CA).	03.12.9	BY, CA, CH, CN, CU, CZ, DB, DK, EB, ES, FI, GB, GE, GH, GM, HR, HU, ID, IL, IS, IP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SL, SK, SL, TJ, TM, TR, TT, UA, UG, US, UZ, VN, YU, ZW, ARIPO patent (GH, GM, KE, LS, MW, SD, SZ, UG, ZW), Burasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CY, DE, DK, ES, FI, FR, GB, GR, CT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CY, DE, DK, ES, FI, FR, GB, GR, CT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CY, DE, DK, ES, FI, FR, GB, CR, CY, DE, DK, ES, FI, FR, GB, GR, CT, CY, DE, DK, ES, FI, FR, GB, GR, CT, CY, DE, DK, ES, FI, FR, GB, GR, CY, CY, CY, DE, DK, ES, FI, FR, GB, GR, CY, CY, CY, CY, CY, CY, CY, CY, CY, CY
(72) Inventors; and (75) Inventors/Applicants (for US only): DOUGLAS, I [CA/CA]; #105-876 W. 16th Avenue, Vancouve Columbia V5Z 1T1 (CA). CAMPBELL, Jenu [CA/US]; Apartment 1, 37 Francesca Street, St MA 02143 (US). COLLINGS, Bruce, A. [CA/C Bainbridge Avenue, Burnaby, British Columbia (CA). (74) Agent: BERESKIN & PARR; 40th floor, 40 King St Toronto, Ontario M5H 3Y2 (CA).	claims and to be republished in the event of the receipt of the amendments. amendments.	

(54) THIE: METHOD OF ANALYZING IONS IN AN APPARATUS INCLUDING A TIME OF FLIGHT MASS SPECTROMETER AND A LINEAR ION TRAP

(57) Abstract

A method of analyzing ions is carried out in a mass spectrometer apparatus comprising an ion source, linear RF quadrupole and a time of flight mass spectrometer. Ions are generated from the ion source and passed into the linear RF quadrupole. To retain ions within the linear RF potentials are quadrupole, applied to either end of it and it is then operated as an ion trap. Ions of interest are selected in the linear RF quadrupole and unwanted ions are caused to be ejected. Selected ions are then excited and caused to collide with a neutral gas, to cause collision induced dissociation thereof, thereby forming fragment ions for analysis in the time of flight mass spectrometer. The potential of one end of the linear RP quadrupole is then adjusted to pass selected and fragment ions through to the

mass spectrometer.

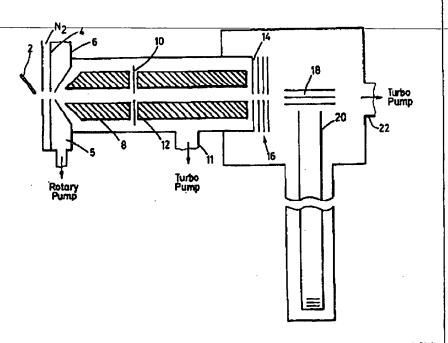

time of flight mass spectrometer. This enables a spectrum of the selected and the fragment ions to be obtained from the time of flight

POR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

Athenia	R.S	Snain	1.9	Lesotho	SI	Slovenia	
• • • • • • • • • • • • • • • • • • • •					5K	Slovakia	
					SN	Senega)	
		• • • • • • • • • • • • • • • • • • • •					ı
			MC				
			MD		TG		-
			MG	•	TJ		i
	GN	Guinea	MK		TM	Torkmenistan	- 1
	GR	Greece			TR	Turkey	ı
Boteuria	80		ML	Mali	TT	Trinidad and Tobago	
Benin	IR	Ireland	MN	Mongolia	UA	Ukraine	ı
Brazil		Israel	MR	Mauritania		Uganda	
Belarus	18	Iceland	MW	Malawi	US	United States of America	- 1
Canada	IT	Italy	MX	Mexico	UZ	Uzbekistan	1
Central African Republic	JP	Japan	NE	Niger	VN	Viet Nam	- 1
Congo	KE		NL	Netherlands	¥U	Yugoslavia	
boston bost w 2	KC		NO	Norway	2W	Zimbabwe	
Côte d'Ivoire	KP		NZ	New Zealand			1
Cameroon		Republic of Korea	PL	Poland			J
Chine	KR	Republic of Korea	PT	Portugal			ı
Cuba	ΚZ	Kezakstan	RO	Romania			- 1
Czech Republic	LC	Saint Lucia	RU	Russian Federation			- 1
Gemany	LJ	Liechtenstein	SD-	Sudan			
Demnark	LK	Szi Lenka	SE	Sweden			1
Estonia	LR	Liberia	SG	Singapore			
							- 1
	Bosnia and Herzogovina Barbados Belgium Barbtha Faso Bolgaria Benin Brazii Belarus Canada Central African Republic Congo Switzerland Côte d'Ivoire Cameroon China Cuba Cuba Centraly	Amenia FI Amenia FR Amenia FR Amenia FR Amenia FR Amenia GA Amenia GB Borola and Herzegovine GE Barbados GH Belgium GN Burtha Faso GR Bolgaria HU Benin IR Brazil IIL Bolana IS Canada IY Contral African Republic JP Congo KE Switzerland KC Chte d'Ivoire KP Cameroon China KR Cuba KZ Czech Republic LC Germany LI Denmant LI Denmant LI Denmant LI Denmant LI LI Denmant LI LI CAR Cuba CL Commant LI Denmant LI LI Denmant LI LI Denmant LI L	Armenia PI Phland Assuria PR Prance Assurila GA Gabon Azorbaijim GB Ualted Kingdom Barbados GH Ghana Belgium GN Gulnea Burkha Paso GR Greece Botgaria BU Hungary Benin IR Irchand Brazil IL Israel Bolana IS Iceland Canada IY Italy Central African Republic JP Japan Congo KS Kenya Swizzerland KC Kyrgyzatan Che d'Ivoire KP Deneceraic People's Cameroon Chiaa KR Republic of Korea Cuba Czech Rapublic LC Saiat Lucia Cermany LI Lockenstoin Denmank LK St Larka	Armenia FI Finland LT Assurals FR Prance LU Ansuralis GA Gabon LV Armenia GB United Kingdom MC Bosnia and Herzogovina GE Georgia MD Barbados GH Ghana MG Belgium GN Gulnea MK Burthas Faso GR Greece Botgaria BU Hungary ML Benin IR Irchard MN Brazil IL Israel MR Bolarus IS keland MW Canada IV Italy MX Canada IV Italy MX Canada IV Italy MX Canada KR Kenya NL Switzerland KG Kyrgyzatas NO Che d'Ivoire KP Democratic People's NZ Cameroon Republic of Korea PL Cuba KZ Kazakstas RO Cecch Republic LC Saiat Lucin RU Germany LI Liechsenstein SD Demnark LK SI Lerks SR	Armenia FI Pinland LT Lithuania Assuria FR Prance LU Lunembourg Australia GA Gabon LV Latvia Aparbaijim GB Unleed Kingdom MC Monaco Bosola and Herzegovina GE Georgia MD Republic of Moldova Barbados GH Chana MG Madagascar Belgium GN Gulnea MK The former Yugostav Berkha Faso GR Greece Republic of Macedonia Botgaria BU Hungary ML Mall Benin IR Irchnd MN Mongolia Brazii IL Israel MR Maurikania Belana IS Icchnd MW Makwi Canada IY Italy MX Mexico Congo KE Kenya NL Netherlands Switzerland KC Kyrgyzstan NO Norway Che d'Ivoire KP Denocraic People's NZ New Zealand Cuba KZ Kazakstan RO Romania Cuba KZ Kazakstan RO Romania Commany LI Lichenestoin SE Sudan Dennark LK Si Larka SE Sweden	Armenia FI Pinland LT Lithuania SK Assuria FR Prance LU Lucembourg SN Australia GA Gabon LV Latvia SZ Aparbatijam GB Ualeed Kingdom MC Monaco TD Bosola and Herzegovina GE Georgia MD Republic of Moldova TG Barbados GH Chana MG Madagascar TJ Belgium GN Gulnea MK The former Yugostav TM Burkhas Faso GR Greece Republic of Macedonia TR Bobgaria BU Hungary ML Mail TT Benin IR Irchnd MN Mongolia UA Brazzii IL Israel MR Mauritania UG Bolarus IS Iccland MW Makawi US Canada IY Italy MX Mexico UZ Central African Republic JP Japan NE Niger VN Congo KE Kenya NL Metheriands YU Switzerland KG Kyrgyzstan NO Norway ZW Cète d'Ivoire KP Deneocratic People's NZ Cuba KZ Kazakstan RO Romania Cuba KZ Kazakstan RO Romania Cuch Republic LC Saiat Lucis RU Rustan Federation Dennark LK Sti Lenka SE Sweden	Armenia FI Finland LT Lithunia 5K Slovatia Assuria FR Prace LU Lunembourg SN Senegal Assurials GA Gaboo LV Lavis SZ Swaziland Armenia GB Unlood Kingdom MC Momaco TD Chad Bonia and Herzegovins GE Georgia MD Republic of Moldova TG Togo Barbados GH Ghana MG Madagascar TJ Tajlikistan Belgium GN Gulosa MK The former Yugoslav TM Torkmenistan Burkha Faso GR Greece Republic of Macdonis TR Turkey Bolgaria BU Hungary ML Mali TT Trinidad and Tobingo Benin IR Ireland MN Mongolia UA Ukrains Brazil IL Israel MR Maurikania UG Uganda Bolana IS Iceland MW Makavi US United States of America Canada IY Italy MX Mexico UZ Uzbekistan Congo KS Kenya NL Netherlands YU Yugoslavia Switzenland KG Kyrgyzstan NO Norway ZW Zimbabwe Chee d'Ivoire KP Deneocratic People's NZ New Zealand Cuba KZ Kazakstan RO Romania Cuba KZ Kazakstan RO Romania Coemany LI Lichenstein SE Sweden

Title: METHOD OF ANALYZING IONS IN AN APPARATUS INCLUDING A TIME OF FLIGHT MASS SPECTROMETER AND A LINEAR ION TRAP

FIELD OF THE INVENTION

This invention relates to mass spectrometry, and more particularly is concerned with a method of analyzing ions in an apparatus including a time of flight mass spectrometer.

BACKGROUND OF THE INVENTION

In the field of mass spectrometry, a wide variety of different spectrometers have been developed, and combinations of different 10 spectrometer elements. One well-known type of spectrometer is a quadrupole mass spectrometer, and it is known to provide devices with two or more mass analyzing quadrupole stages to provide MS/MS capabilities. It is also known to combine a quadrupole stage with a time of flight mass spectrometer (TOF MS), as detailed below. A TOF MS has the advantages of high scan speed, unlimited mass range and, if a reflectron is used, a resolution of 10,000 or more. However, TOF MS does not normally provide MS/MS capabilities. 3D ion trap mass spectrometers can perform MS/MS analysis in a comparatively simple device, but generally are operated at a lower resolution than a reflectron TOF MS. Higher resolution can be 20 achieved with an ion trap mass spectrometer, but only with very low scan speeds. As well, it is difficult to inject ions from an external source into a 3D ion trap and the mass range is limited.

One proposal by one of the inventors of the present invention is disclosed in U.S. Patent 5,179,278. This patent describes the use of an RF multipole ion guide as an interface between an ion source and an ion trap. The intention is to improve the duty cycle of the ion trap mass spectrometer. However, there is no specific teaching of using the multipole device itself in a full trapping mode to provide MS/MS capabilities. Rather,

the teaching is that the trapping is obtained by applying selected electric potentials at the ends of the space in the multipole device, to cause ions to be reflected from a second, outlet end to a first, inlet end and then back again towards the second, outlet end. This retains the ions in the space for a longer period of time, this being long enough for analysis to be carried out on a previous ion sample fed to the ion trap spectrometer.

U.S. Patent 5,652,427 describes the use of an RF multipole that includes a number of separate stages with different degrees of vacuum. However, the arrangement is intended simply to transfer ions from a high pressure source to a mass spectrometer. There is no teaching of trapping of ions in any multipole stage nor any teaching of MS/MS capabilities through resonant excitation and ejection and the like.

U.S. Patent 5,420,425 (Bier et al. and assigned to Finnigan Corporation) is concerned with an ion trap mass spectrometer, for analyzing 15 ions. It has electrodes shaped to promote an enlarged ion occupied A quadrupole field is provided to trap ions within a volume. predetermined range of mass-to-charge ratios, and the field is then changed, so that trapped ions with specific masses become unstable and leave the trapping chamber in a direction orthogonal to the central axis of the chamber. The ions leaving the spectrometer are detected, to provide a signal indicative of their mass-to-charge ratios. The patent does teach a method of first introducing ions within a pre-determined range of mass-tocharge ratios into the chamber and subsequently changing the field to select just some of the ions for further manipulation. The quadrupole field is then adjusted so as to be capable of trapping product ions of the remaining ions. The remaining ions are then dissociated or reacted with a neutral gas to form those product ions. The quadrupole field is then changed again, to remove, for detection, ions whose mass-to-charge ratios lie within the desired range. It is noteworthy that the ions are not detected by a Time of

Flight (TOF) instrument. The Finnigan device, since it uses radial ejection, will produce a stream of ions having broad space and velocity distributions. It would be difficult to manage such a beam and introduce it into a TOF MS analyzer.

Other workers have demonstrated interfacing a TOF MS to an electrospray source, using a linear RF quadrupole operated at moderately high pressure (for example Chemushevich et al., presented at the 44th ASMS Conference on Mass Spectrometry and Allied Topics, Portland, Oregon, May 12-16, 1996). Also, other workers have used hexapole and 10 octopole ion guides instead of a quadrupole. In hexapole and octopole fields, ions of different m/z do not normally have well defined frequencies of motion and so resonant excitation or ejection of selected ions is not possible, and this is a significant advantage of using a quadrupole.

It has also been known to use a 3D ion trap as an interface 15 between an ion source and a TOF MS (S. M. Michael et al, Rev. Sci. Instr. 63, 4277-4284, 1992; Purves and Li, J. Microcolumn Separations 7 (6) 603, 1995). The 3D ion trap can be provided with MS/MS capability (Qian and Lubman, Rap. Commun. Mass. Spec. 10, 1079, 1996). Use of a three dimensional ion trap has a number of disadvantages. Firstly, ion injection efficiency is at 20 least ten times less than the efficiency with a two dimensional quadrupole. Secondly, the ion storage volume is less in the three dimensional trap, so that only a relatively small number of ions can be stored, without space charge problems, and for this reason, the concentration dynamic range is limited in a three dimensional trap.

A related approach has been proposed using two separate multipoles and then a TOF mass spectrometer (H. R. Morris et al., Rap. : Commun. Mass. Spec. 10, 889, 1996). Here, selection of ions in a given m/z is carried out conventionally in a first quadrupole mass filter. These are then passed to an RF only hexapole and dissociated by collision with a

neutral gas. Resulting ions then pass through to a TOF MS, to obtain a spectrum of the product ions. In addition, a system with a first mass analyzing quadrupole and a second RF only quadrupole has been described (Chevchenko Rapid Communications in Mass Spectrometry, Vol 11, Page 1015-1024, 1997). Both of these systems are relatively complex and expensive and have a number of stages, which will likely lead to loss of sensitivity.

Finally, a more recent proposal is found in an article entitled "A New Technique for Decomposition of Selected Ions in Molecule Ion Reactor Coupled with Ortho-Time-of-Flight Mass Spectrometry" (A. 10 Dodonov et al. Rapid Communications in Mass Spectrometry, 11, 1649-1656, 1997). This paper shows restricted experimental results carried out on, in effect, pre-selected ions, i.e. the tests were carried out with a single chemical compound. There is no specific teaching of using the apparatus to carry out the selection stage. Two modes of ion dissociation are disclosed. In a first mode, the motion of parent and fragment ions is chosen to be stable, and the RF electric field forces the ions to oscillate around the quadrupole axis. Simultaneously, a DC potential is applied along the axis of the device to accelerate the ions, and the strength of this field controls the collision induced fragmentation of the ions. The quadrupole is filled with gas at around 1 mbar of pressure for this purpose. There is no ion mass to charge ratio selection and all ions present will be accelerated by this field. Fragment ions with mass-to-charge ratios both above and below the m/z of a parent ion can be transmitted to a TOF mass analyzer for analysis. It can be noted that the applied field will also cause acceleration and possible further 25 fragmentation of the fragment ions, although control over the strength of the field can limit this to some extent. Nonetheless, the applied axial field will not discriminate between the different types of ions.

In a second mode, ions are fragmented by confining them in the quadrupole with the RF electric field chosen to have an amplitude

and frequency such that desired ions are close to the limit for stable ion motion with Mathieu's parameter q=0.9. This causes an increase of the parent ion velocity and thus leads to collisional "heating" and fragmentation of the parent ions. Then, only those ions with m/z ratios above that of the precursor or parent ion are stable in the quadrupole and only these ions are transmitted to the detector, i.e. the TOF MS. Ions with an m/z ratio smaller than that of the parent ion are rejected due to the unstable character of their motion. For multiply charged ions this may not be a severe limitation because some fragment ions may have lower charges and hence higher mass to charge ratios, but for singly charged ions no fragments will be detected.

A disadvantage of both modes is that all ions in the quadrupole are excited at the same time and dissociate. If there are two compounds present, they would both fragment and in general it would not be possible to tell which fragments came from which precursor.

SUMMARY OF THE INVENTION

The present inventors have realised that one can obtain the capabilities of a tandem mass spectrometer in a relatively simple device, by combining a linear quadrupole, or other multipole, with a TOF MS. The quadrupole, or other multipole, is operated as an ion trap, and an ion is selected by resonant ejection of ions of other masses or otherwise. The isolated ions are then excited and caused to undergo collision induced dissociation or fragmentation, in the quadrupole or other multipole.

In accordance with the present invention, there is provided a method of analyzing ions in a mass spectrometer apparatus comprising an ion source, a linear RF quadrupole and a time of flight mass spectrometer, the method comprising the steps of:

(1) generating ions from the ion source and passing the

ions into the linear RF quadrupole;

- (2) applying potentials at either end of the linear RF quadrupole and operating the linear RF quadrupole as an ion trap;
- (3) selecting ions of interest in the linear RF quadrupole and ejecting unwanted ions;
 - (4) exciting the ions and causing the ions to collide with a neutral gas, to cause collision induced dissociation thereof, thereby forming fragment ions for analysis in the time of flight mass spectrometer.;
- (5) adjusting the potential of one end of the linear RF quadrupole, to pass the selected and the fragment ions through to the time of flight mass spectrometer; and
 - (6) obtaining a spectrum of the selected and the fragment ions in the time of flight mass spectrometer.

It is preferred to use a quadrupole device, as these 15 inherently have well-defined stability parameters and frequencies of excitation for a particular ion. The x and y motions are separate and either can be excited with good selectivity. However, for some applications, it might be desirable or possible to use other 2D multipole designs, such as a hexapole or octopole. If such an instrument is operated with low RF 20 voltages, the ion motion is approximately harmonic motion with well defined frequencies (as described by Gerlich in Advances in Chemical Physics vol 82 1992, pages 1-176). It should also be noted that the linear RF quadrupole or multipole can comprise a single quadrupole or multipole or alternatively two sets of quadrupole or multipole rods can be provided in tandem.

A further aspect of the present invention is that the linear RF quadrupole can be used to carry out multiple mass spectrometry steps, so as to perform MSn. Thus the method can include, after step (4), an additional step of isolating and exciting one or more fragment ions, in the

linear RF quadrupole, to cause collision induced dissociation of one or more of the fragment ions to form further fragment ions. Further, the method-can include multiple cycles of isolating and exciting one or more of the fragment ions in the linear RF quadrupole, wherein each cycle comprises isolating and exciting at least one or more of the fragment ions formed in the previous cycle to form further fragment ions. In each cycle all the fragment and the selected ions can be excited to cause collision induced dissociation.

The selected ions and/or the fragment ions can be excited

by one of (i) exciting the selected ions by resonance excitation at a particular
secular frequency and (ii) applying a broadband excitation waveform, to

cause collision induced dissociation of the selected ions.

Yet another aspect of the invention provides an apparatus incorporating one or more linear RF quadrupoles, or other multipoles, and a time of flight mass spectrometer, and adapted to carry out the method of the present invention.

BRIEF DESCRIPTION OF THE DRAWING

For a better understanding of the present invention and to show more clearly how it may be carried into effect, reference will now be made, by way of example, to the accompanying drawing, which shows schematically a preferred embodiment of the present invention.

DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENT

With reference to the drawing, an electrospray source is indicated at 2. It will be appreciated that any suitable ion source could be used such as EI (electron ionization), CI (chemical ionization), laser desorption, etc. Ions from the source 2 pass through an orifice 4. A supply of nitrogen gas, which maybe heated, is provided as indicated, to promote

vaporization of solvent. The ions then pass into a chamber 5 provided with a connection to a rotary pump, to maintain a desired low pressure. A skimmer 6 then provides an orifice through which desired ions pass into a first RF quadrupole 8. In known manner, this includes a quadrupole rod set provided with usual connection for supply of RF and DC voltages. The quadrupole 8 is operated in the RF only mode, to transmit ions of a broad range of mass to charge ratios. For simplicity, details of electrical connections, and electrical supplies are omitted.

An entrance lens 10 separates the first quadrupole 8 from a second RF quadrupole 12, but note that the lens 10 does not separate two chambers as the two quadrupoles 8, 12 are essentially in a single chamber, although two chambers at different pressures could be used. The second quadrupole 12 is also operated in the RF only mode. As indicated at 11, a connection is provided to a turbo pump, for maintaining a pressure of, for example, around 1-10 millitorr. As indicated, the first quadrupole 8 is shorter than the second quadrupole 12. For example, the first quadrupole 8 can have a length of 5 cm and the second quadrupole 12 can have a length of 20 cm, i.e. the quadrupoles need not be of the same length.

At the exit from the second quadrupole 12, there is an exit aperture 14 and then a series of additional lens or electrodes 16, for controlling the ion beam and ensuring that it passes into a source region 18 of a time of flight mass spectrometer (TOF MS) 20.

Here, the time of flight mass spectrometer 20 is shown orthogonal to the axis of the quadrupoles 8, 12. It will be appreciated that the TOF MS 20 could equally be axially arranged relative to the quadrupoles 8, 12. In known manner, a connection 22 is provided, to enable the TOF MS to be pumped down to the desired level of vacuum.

In use, if potentials of the entrances and exits of the quadrupoles 8, 12 are set to continuously transmit ions, then conventional

operation of the TOF MS 20 gives a mass spectrum of ions from the source 2. In known manner, the electrodes at the source region 18 of the TOF MS 20 are activated to collect and provide pulses of ions travelling through the TOF MS, whose time of flight is measured, to give a spectrum for those 5 ions.

Now, in accordance with the present invention, stopping potentials can be applied at the entrance and exit of either one or both of the first and second quadrupoles 8, 12. This serves to trap ions in the respective quadrupole. Then, unwanted trapped ions can be ejected by resonant 10 excitation at the secular frequencies of the ions. Also, ions of a single m/z value can be trapped and isolated by ejection of all other ions with a filtered noise field or SWIFT waveform as is known, this essentially being a noise waveform with a notch or gap at the frequency corresponding to the secular frequency of the ion of interest. The isolated ions can then be excited and 15 dissociated by collision with a neutral gas. There are a number of ways of causing excitation, collision and fragmentation. Then, the resultant ion fragments can be driven or transferred into the TOF MS 20 by lowering the trapping voltage on the electrode 14 between the second quadrupole 12 and the TOF MS 20. The ions enter the source region of the TOF MS 20 and the mass spectrum of fragment ions can be obtained.

Ions can be permitted to enter the source region 18 of the TOF MS with near thermal energies that they may have from the second quadrupole 12. Alternatively, they can be accelerated towards the source region 18 by setting up a suitable axial field in the second RF quadrupole 12 25 (as described by B.Thomson et al., at the 44th ASMS conference on mass spectrometry and allied topics, May 12-16th, Portland, Oregon, 1996). Thermal ions typically take of the order of tens of milliseconds to transfer to the source region of the TOF MS and accelerating the ions into the TOF MS 20 has the advantage of reducing the transfer time down to ca. 1 ms.

Thus, for example, in ICP-MS, intense Ar* ions can cause difficulties and can effectively paralyze the detector. To overcome this, a sample of ions from the source 2 can be passed into the first quadrupole 8. There, potentials can be applied to the skimmer 6 and lens 10, and a field applied at the resonance frequency of the Ar* ions to effectively eject them from the ion sample. Then, the voltage on the lens 10 can be adjusted to cause the ion sample to pass into the second quadrupole 12. There, a filtered noise field or SWIFT waveform could be provided to further isolate an ion of interest. The voltage at the aperture 14 and lens 16 would then be adjusted to cause the desired ion to pass into the TOF MS 20.

WE CLAIM:

5

- A method of analyzing ions in a mass spectrometer apparatus comprising an ion source, a linear RF quadrupole and a time of flight mass spectrometer, the method comprising the steps of: generating ions from the ion source and passing the ions (1) into the linear RF quadrupole; applying potentials at either end of the linear RF quadrupole and operating the linear RF quadrupole as an ion trap; selecting ions of interest in the linear RF quadrupole (3) 10 and ejecting unwanted ions; exciting the selected ions and causing the ions to collide with a neutral gas, to cause collision induced dissociation thereof, thereby forming fragment ions for analysis in the time of flight mass spectrometer.; adjusting the potential of one end of the linear RF 15 quadrupole, to pass the selected and the fragment ions through to the time of flight mass spectrometer; and obtaining a spectrum of the selected and the fragment (6)
- A method as claimed in claim 1, wherein step (3) comprises one of (i) 2. 20 ejecting ions having a predetermined secular frequency by excitation at said secular frequency, and (ii) applying a filtered noise field to eject ions other than a desired ion having a single m/z value.

ions in the time of flight mass spectrometer.

A method as claimed in claim 1 or 2, which includes, after step (4), an additional step of exciting one or more of the fragment ions, in the linear RF quadrupole, to cause collision induced dissociation of one or more of the fragment ions to form further fragment ions.

- 4. A method as claimed in claim 3, which includes multiple cycles of exciting one or more of the fragment ions in the linear RF quadrupole, wherein each cycle comprises exciting one or more of the fragment ions formed in the previous cycle to form further fragment ions.
- 5 5. A method as claimed in claim 4, wherein in each cycle all the fragment and the selected ions are excited to cause collision induced dissociation.
- 6. A method as claimed in claim 1, which includes exciting the selected ions by one of (i) exciting the selected ions by resonance excitation at a particular secular frequency and (ii) applying a broadband excitation waveform, to cause collision induced dissociation of the selected ions.
- 7. A method as claimed in claim 5, which includes exciting the selected ions and one or more of the fragment ions by one of (i) exciting the selected and fragment ions by resonance excitation at a particular secular frequency 15—and (ii) applying a broadband excitation waveform, to cause collision induced dissociation of the selected and fragment ions.
 - 8. A method of analyzing ions in a mass spectrometer apparatus comprising an ion source, a linear RF multipole and a time of flight mass spectrometer, the method comprising the steps of:
- 20 (1) generating ions from the ion source and passing the ions into the linear RF multipole;
 - (2) applying potentials at either end of the linear RF multipole and operating the linear RF multipole as an ion trap;
 - (3) selecting ions of interest in the linear RF multipole and


ejecting unwanted ions;

- (4) exciting the selected ions and causing the ions to collide with a neutral gas, to cause collision induced dissociation thereof, thereby forming fragment ions for analysis in the time of flight mass spectrometer.;
- (5) adjusting the potential of one end of the linear RF multipole, to pass the selected and the fragment ions through to the time of flight mass spectrometer; and
- (6) obtaining a spectrum of the selected and the fragment ions in the time of flight mass spectrometer.
- 9. A method as claimed in claim 8, wherein step (3) comprises one of (i) ejecting ions having a predetermined secular frequency by excitation at said secular frequency, and (ii) applying a filtered noise field to eject ions other than a desired ion having a single m/z value.
- 15 10. A method as claimed in claim 8 or 9, which includes, after step (4), an additional step of exciting one or more of the fragment ions, in the linear RF multipole, to cause collision induced dissociation of one or more of the fragment ions to form further fragment ions.
- 11. A method as claimed in claim 10, which includes multiple cycles of exciting one or more of the fragment ions in the linear RF multipole, wherein each cycle comprises exciting one or more of the fragment ions formed in the previous cycle to form further fragment ions.
- 12. A method as claimed in claim 11, wherein in each cycle all the fragment and the selected ions are excited to cause collision induced dissociation.

- 13. A method as claimed in claim 8, which includes exciting the selected ions by one of (i) exciting the selected ions by resonance excitation at a particular secular frequency and (ii) applying a broadband excitation waveform, to cause collision induced dissociation of the selected ions.
- 5 14. A method as claimed in claim 12, which includes exciting the selected ions and one or more of the fragment ions by one of (i) exciting the selected and fragment ions by resonance excitation at a particular secular frequency and (ii) applying a broadband excitation waveform, to cause collision induced dissociation of the selected and fragment ions.
- 10 15. A method as claimed in claim 8, which comprise carrying out the method using one of a hexapole and an octopole as the linear RF multipole.

PCT/CA98/01106

WO 99/30350 PCT/CA9

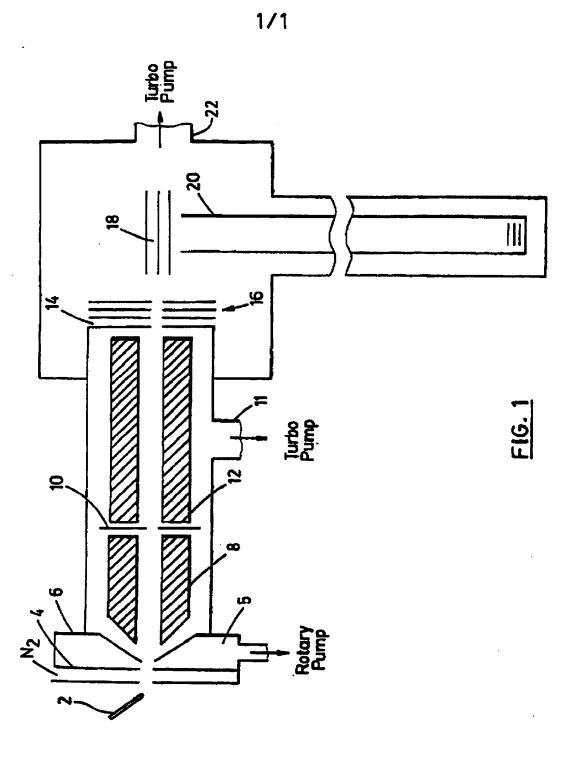

SUBSTITUTE SHEET (RULE 26)



INTERNATIONAL SEARCH REPORT

ti rational Application No PCT/CA 98/01106

According	to international Patent Classification (IPC) or to both national classifica-	illon end (PC	
B. FIELD	S SEARCHED		
Minimum IPC 6	documentation searched (classification system followed by classification HO1J	n symbola)	
Documen	ation secreted other than minimum documentation to the extent that se	uch documents are included in the field	de esarched
Electronic	data base consulted during the international search (name of data bas	e and, where practical, search terms	used)
C. DOCU	MENTS CONSIDERED TO BE RELEVANT		
Category	Citation of document, with indication, where appropriate, of the rele	vard passages	Relevant to claim No.
Y	US 5 689 111 A (DRESCH THOMAS ET 18 November 1997		1-15
	see column 6, line 6 - column 6, figure 5	(
Y	US 5 576 540 A (JOLLIFFE CHARLES 19 November 1996		1-15
	see column 8, line 62 - column 9, figure 1		
	see column 10, line 32 - column 10		
	see column 10, line 53 - column 10 67 see column 11, line 21 - column 11		
	26		
	-/	/ 	
	·		
X Furt	her documents are listed in the continuation of box C.	X Patent family members are list	ed in annex.
Special categories of cited documents: "A" document defining the general state of the art which is not considered to be of particular relevance. "E" earlier document out published on or after the international		later document published after the it or priority date and not in conflict wi cited to understand the principle or invention	theory underlying the
"L" docume which dtation	of which may throw doubte on priority claim(s) or	 document of particular relevance; the cannot be considered novel or can involve an inventive step when the document of particular relevance; the cannot be considered to involve an 	not be considered to document is taken alone e claimed invention
P docume	office means document is combined with one or n ments, such combined with one or n ments, such combined with one or n ments, such combined by in the st. "å" document member of the same pater		more other such docu- rious to a person sisted
Date of the	actual completion of the international search	Date of mailing of the international s	•
	5 April 1999	18/05/1999	
Name and n	alling address of the ISA European Patent Office, P.B. 5816 Patentiaan 2	Authorized officer	

INTERNATIONAL SEARCH REPORT

			PCT/CA 00/01106	
	C.(Coatinu	untion) DOCUMENTS CONSIDERED TO BE RELEVANT	PCT/CA 98/01106	
	Category *	Chatton of document, with Indication, where appropriate, of the relevant passages	Relevant to claim No.	-
	A -	WO 97 07530 A (MDS HEALTH GROUP LTD; THOMSON BRUCE A; JOLLIFFE CHARLES L) 27 February 1997 see page 26, line 10 - page 26, line 14	1-15	
	A	GB 2 300 751 A (BRUKER FRANZEN ANALYTIK GMBH) 13 November 1996 see page 9, paragraph 1-3	1-15	
	A	EP 0 529 885 A (MDS HEALTH GROUP LTD) 3 March 1993 cited in the application see the whole document	1,8	
	A	US 5 420 425 A (BIER MARK E ET AL) 30 May 1995 cited in the application see the whole document	1,8	
	A	CHIEN B M ET AL: "ANALYSIS OF THE FRAGMENTS FROM COLLISION-INDUCED DISSOCIATION OF ELECTROSPRAY-PRODUCED PEPTIDE IONS USING A QUADRUPOLE ION TRAP STORAGE/REFLECTRON TIME-OF-FLIGHT MASS SPECTROMETER" ANALYTICAL CHEMISTRY, vol. 66, no. 10, 15 May 1994, pages 1630-1636, XP000449869	1,8	
	A	SHEVCHENKO A. ET AL: "Rapid 'de Novo' Peptide Sequencing by a Combination of Nanoelectrospray, Isotopic Labeling and a Quadrupole/Time-of-Flight Mass Spectrometer" RAPID COMMUNICATIONS IN MASS SPECTROMETRY, vol. 11, 1997, pages 1015-1024, XP002101143 cited in the application see figure 1	1,8	,
4				
Fı	im PCT/ISA/210 (oor	rifinuation of second ether() (July 1992)		

page 2 of 2

Ir. attonal Application No

	INTE	ERNATIONAL SEAR	CH REPORT		
	-	information on patent family		in ational Application No	
	Patent document clied in search report	Publication date	Patent famil		
	US 5689111	A 18-11-1997	WO 9807:	597 A 06-03-1998 177 A 19-02-1998 178 A 19-02-1998	
	US 5576540	A 19-11-1996	NONE	19-02-1998	
	WO 9707530	A 27-02-1997	AU 66532 CA 22290 EP 08438 US 58473	070 A 27-02-1997 387 A 27-05-1998	
	GB 2300751	A 13-11-1996	DE 195175 US 57082	507 C 08-08-1996	
	EP 0529885	A 03-03-1993	US 51792	278 A 12-01-1993 28 A,C 24-02-1993 23 D 18-01-1995	:
	US 5420425	A 30-05-1995	CA 21483 DE 695026 DE 695026 EP 06846 JP 26580 JP 732632	31 A 28-11-1995 62 D 02-07-1998 62 T 24-12-1998 28 A 29-11-1995 12 B 30-09-1997	
·				12 12 1373	
•					
		_			
Fa	em PCT/ISA/210 (patent family ennex) (July 15	702)			

This Page is Inserted by IFW Indexing and Scanning Operations and is not part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

BLACK BORDERS

IMAGE CUT OFF AT TOP, BOTTOM OR SIDES

FADED TEXT OR DRAWING

BLURRED OR ILLEGIBLE TEXT OR DRAWING

SKEWED/SLANTED IMAGES

COLOR OR BLACK AND WHITE PHOTOGRAPHS

GRAY SCALE DOCUMENTS

LINES OR MARKS ON ORIGINAL DOCUMENT

REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.